

Super-Snap SuperBuff

CLICK / SCAN



FOR MORE INFO

IMPREGNATED FELT DISKS FOR SUPER POLISHING DIRECT RESTORATIONS

Designed for clinicians who prefer the convenience of achieving a high gloss super polish without the mess of polishing pastes

- o Single use, snap-on, aluminum oxide impregnated felt disks
- o Fast, convenient and comfortable for the patient
- o Ideal to polish natural teeth after scaling

Super-Snap SuperBuff felt disks are designed for single use and can be used with the same CA mandrel as the Super-Snap disk system. Proven to provide a high lustre super polish similar to natural enamel for enhanced



Available in:


- **SuperBuff Set (REF 0535)**
25x unmounted polishers and 2 mandrels (CA shank)
- **Refills**
50x unmounted polishers [REF 0536]
6x CA Mandrels [REF 0439]

Super-Snap
SuperBuff



REF 0536

Super-Snap
Mandrel CA

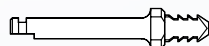

REF 0439


CLINICAL APPLICATIONS

- High gloss polishing of all tooth coloured direct restorative materials
- Polishing after routine scaling

FEATURES

- ✓ **Mess -free and convenient to use**
 - Snap-on impregnated felt disk
 - Soak in water for 30sec before use
 - Same Mandrel as Super-Snap system
- ✓ **Fast, safe and predictable super polishing**
 - Thick disk design enables super polishing of occlusal surfaces efficiently
 - Soft, safe, elastic silicone snap-on for easy placement on mandrels
 - No metal centre to completely eliminate the risk of gouging restoration
- ✓ **Dry polishing process**
 - Low speed (max 15.000 rpm) with uniform and feather light pressure
- ✓ **Disposable disks** for enhanced infection control

INSTRUCTION FOR USE

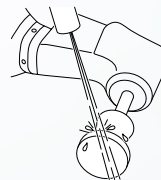

Wet the surface to be polished and moisten the SuperBuff disk with water

RELATED PRODUCTS

- OneGloss 2-in-1 Smart Polisher
- Super-Snap X-Treme Disk System
- Beautifil II / Beautifil II LS/ Beautifil Injectable X



Official Partner

